---

# **Instructions for Use**

## **Diagnostic Suite 2.1.0**

---

|           |                                   |
|-----------|-----------------------------------|
| <b>EN</b> | <b>Instructions for Use</b>       |
| <b>DE</b> | <b>Gebrauchsanweisung</b>         |
| <b>SV</b> | <b>Bruksanvisning</b>             |
| <b>FI</b> | <b>Käyttöohje</b>                 |
| <b>NO</b> | <b>Bruksanvisning</b>             |
| <b>DA</b> | <b>Brugsanvisning</b>             |
| <b>NL</b> | <b>Gebruiksaanwijzing</b>         |
| <b>FR</b> | <b>Instructions d'utilisation</b> |
| <b>IT</b> | <b>Istruzioni per l'uso</b>       |
| <b>ES</b> | <b>Manual de instrucciones</b>    |
| <b>PT</b> | <b>Instruções de Utilização</b>   |
| <b>PL</b> | <b>Instrukcja obsługi</b>         |
| <b>CS</b> | <b>Návod k použití</b>            |
| <b>HU</b> | <b>Használati útmutató</b>        |
| <b>RO</b> | <b>Instrucțiuni de utilizare</b>  |
| <b>EL</b> | <b>Οδηγίες χρήσης</b>             |
| <b>ET</b> | <b>Kasutusjuhend</b>              |
| <b>RU</b> | <b>Руководство по применению</b>  |

# **Instructions for Use - EN**

**Interacoustics Diagnostic Suite 2.1.0**

# Table of contents

|       |  |    |
|-------|--|----|
| 1     | Introduction .....   | 1  |
| 2     | System Requirements.....   | 1  |
| 3     | Installation and System Setup .....                                  | 2  |
| 4     | Starting Diagnostic Suite.....                                       | 3  |
| 4.1   | Start from Noah .....  | 3  |
| 4.2   | Start from OtoAccess™ .....  | 4  |
| 4.3   | Start as Standalone .....  | 4  |
| 5     | Diagnostic Suite – instrument configuration .....                    | 5  |
| 5.1   | Instrument Setup .....   | 5  |
|       | 5.2 Transferring Measurement Data.....                               | 7  |
| 6     | Diagnostic Suite with the New AD226/AD229/AD629/ AC40/AT235(h) ..... | 8  |
| 6.1   | Instrument setup .....   | 8  |
| 6.2   | General setup .....  | 9  |
| 6.3   | Suite Settings.....  | 10 |
| 6.4   | Protocol setup.....  | 11 |
| 6.5   | SYNC Mode.....   | 12 |
| 6.4.1 | One click data transfer (Hybrid Mode disabled) .....                 | 12 |
| 6.4.2 | The Sync Tab .....   | 12 |
| 6.5.1 | 6.4.3 Client Upload .....  | 12 |
| 6.4.4 | Session download .....   | 13 |
| 6.6   | Quick Start – Transferring and Saving Tympanometry Data.....         | 14 |
| 6.7   | Hybrid (Online/PC-Controlled) Mode .....                             | 15 |
| 6.7.1 | Using the Tone Screen.....   | 15 |
| 6.7.2 | Using the Speech Screen.....   | 22 |
| 6.7.3 | Speech Audiometry in Graph Mode .....                                | 24 |
| 6.7.4 | Speech Audiometry in Table Mode .....                                | 25 |
| 6.7.5 | PC Keyboard Shortcuts Manager.....                                   | 28 |
| 7     | Using the Print Wizard .....   | 29 |

## 1 Introduction

This manual describes the Diagnostic Suite software used for transferring audiometric and tympanometry data from the Interacoustics standalone audiometers to the PC. The Diagnostic Suite allows users to display, store and print audiometric data.

## 2 System Requirements

### General computer requirements

- 2 GHz Intel Core 2 Duo CPU
- 2GB Ram
- 1.5 GB available disk space

### Display requirements

- 1024x768 resolution
- Hardware accelerated DirectX/Direct3D graphics card.

### Software requirements

- Windows® XP (SP3), Windows® Vista SP3 (x86),
- Windows® 7 SP1 (x86 and x64)
- Windows 8 / 8.1 (x86 and x64)

- Database support:
  - OtoAccess™
  - Noah 3.7 or Noah 4 (from HIMSA)
  - AuditBase System 4
  - Mirage
  - Practice Navigator
  - Power Office
  - AkuWin
  - He@r-O
  - Entendre 2
  - Quasar
- Instrument to PC connection support:
  - AS608e, AD226 hybrid, new AD629(AD229)-b/e hybrid, new AC40 hybrid, New AT235(h) uses a direct USB connection (it has build-in USB)
  - Old AT235/AA222 has a build-in UCA40 with a USB connector on the rear panel, follow the UCA40 installation procedure.
  - Old AD226, old AD229b/e, Old AC40, use the UCA40 USB-serial converter/adapter.
- Test measurement data support:
  - Audiometric data: Air, Bone, Speech
  - Tympanometry data: Tympanogram, Acoustic Reflex, Eustachian tube function test for the non-perforated eardrum (ETF1) and perforated eardrum (ETF2)

### 3 Installation and System Setup

The Diagnostic Suite can be installed with either OtoAccess<sup>TM</sup> or Noah or run as a standalone application.

To use the software in conjunction with a database (e.g. Noah3.7, Noah4 or OtoAccess<sup>TM</sup>), make sure the database is installed prior to the Diagnostics Suite installation. Follow the manufacturer's installation instructions provided to install the relevant database.

Note that if you are using AuditBase System 4 you must ensure to launch this office system prior to this installation of the Diagnostics Suite.

#### Installation on various Windows® versions

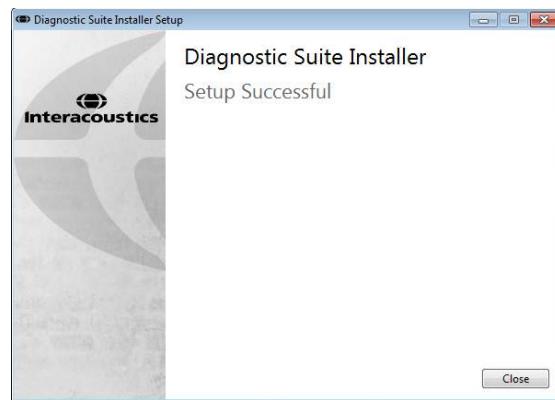
Installation on Windows® XP (SP2 or later), Windows Vista and Windows® 7 (32 and 64 bit) and Windows® 8 systems are supported.

#### Software Installation on Windows® 7

Insert the installation DVD and follow the steps below to install the Diagnostics Suite software. If the installation procedure does not start automatically, click "Start", then go to "My Computer" and double click the DVD/CD-RW drive to view the contents of the installation DVD. Double click the "setup.exe" file to initiate the installation.

Welcome Dialog: Press "Install"

Installation is complete. Press "Close".

Then proceed to system setup described in the next chapter.

Once the Diagnostic suite has been installed, it can be launched from either Noah or OtoAccess<sup>TM</sup> depending on the database you are using.

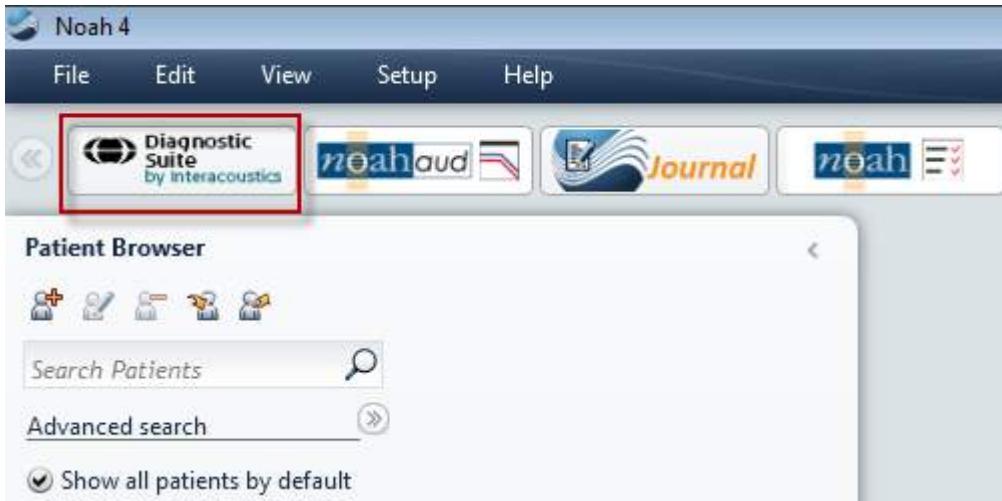
## 4 Starting Diagnostic Suite

The Diagnostic Suite can be run in the following 3 ways:

- 1) From Noah 3 / 4
- 2) From OtoAccess™
- 3) Standalone (without a database)

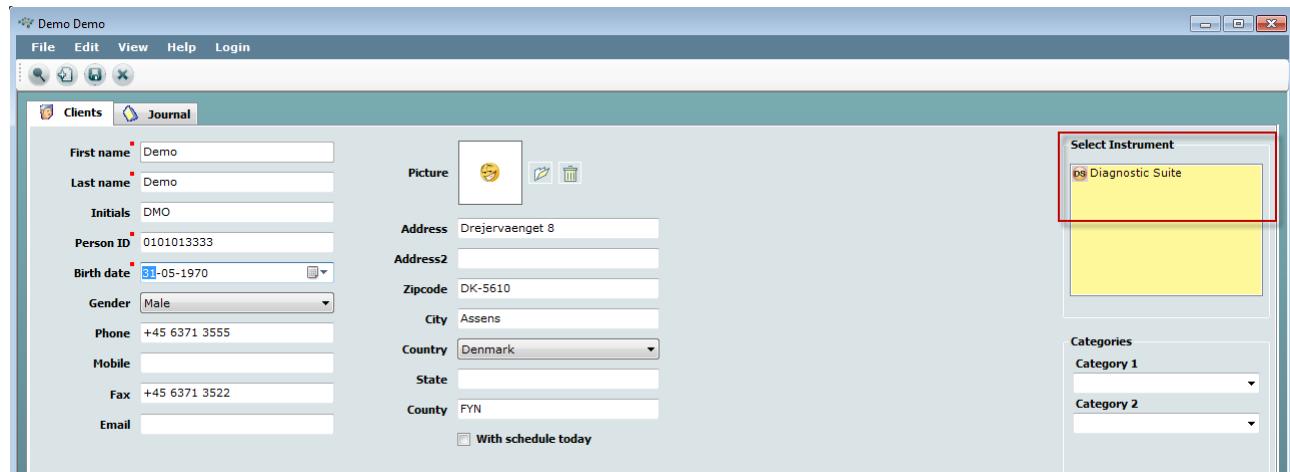
### 4.1 Start from Noah

From Noah activate the Module Selection dialog and select the new Diagnostic Suite icon (below, right side):

## 4.2 Start from OtoAccess™

From OtoAccess™ select the new Diagnostic Suite icon from the “Select Instrument” group box in the top right side:

## 4.3 Start as Standalone

The Diagnostic Suite application may also be run in a standalone mode, without the use of database. In this case the Suite will not handle a patient, only a list of sessions stored in a XML file located under:

*\Documents and Settings\All Users\Application Data\Interacoustics\Diagnostic Suite\Data\Standalone.xml*

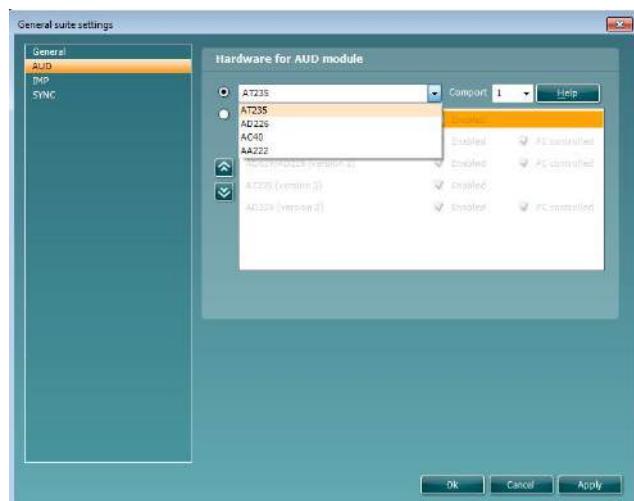
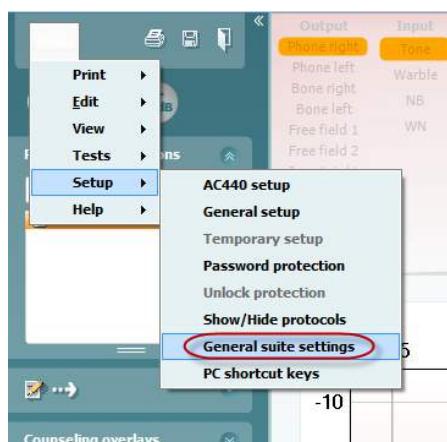
## 5 Diagnostic Suite – instrument configuration

The first section describes how to transfer audiometric data from the instruments that do not support on-line/PC-operated mode: AS608e, AD226, AD229b/e, AC40 (version 1), MT10(v1), New MT10(v2), old AT235 and AA222.

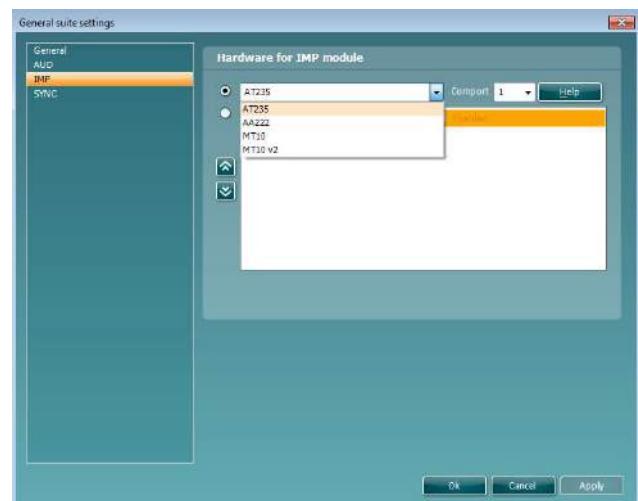
The new AD629/AC40/AD226/AT235(h) supports hybrid mode (on-line / PC-controlled modes) and patient/session transfer which is described in the next chapter.

### 5.1 Instrument Setup

The instrument setup is done in the Diagnostic Suite under **Menu | Setup | General suite settings** under the **AUD / IMP tab** section:

Hardware selection for AUD-module

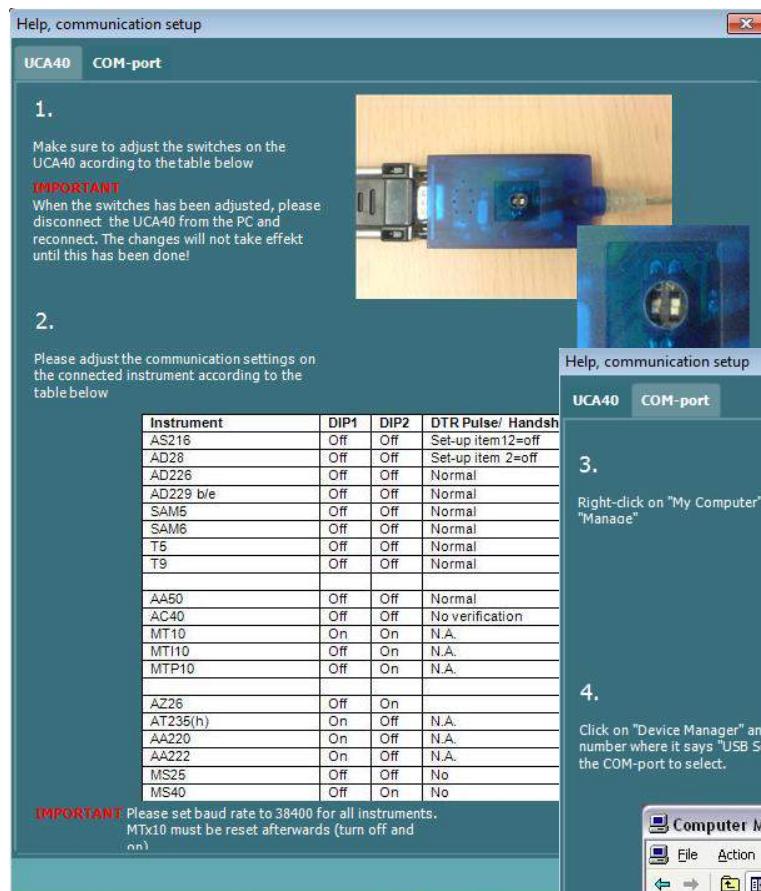
Hardware selection for IMP-module

**Important:** Please be sure not to select the “AD226 (version 2)”, “AD629 (version 2)” or “AC40 (version 2)” because they refer to the new generation USB-based audiometers. Please refer to the next chapter.

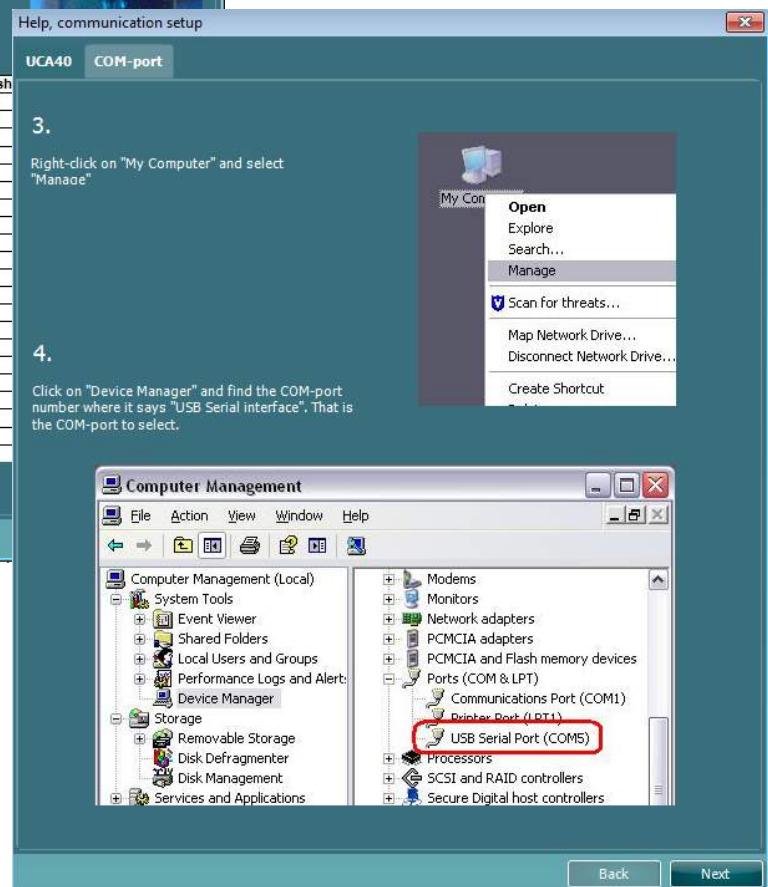
Select from the **Instrument** drop down menu to which instrument you are connected. Then select from the **Comport** drop down menu through which comport the instrument connects to your pc. Note that although your instrument might be connected through a USB connection, you still need to select through which comport that USB connection is supported. To find the appropriate comport, right click on My Computer (on your desktop or in Widows Explorer) and select “Manage”. Click on “Device Manager” and find the comport number where it says “USB Serial Port”. In general the lowest available comport mentioned here is used.

Press **OK** to save your settings and close the General suite settings window.

Press **Help** for further instructions on how to find the appropriate comport and also for the instruments which are connected through an UCA40 to find the appropriate setup.

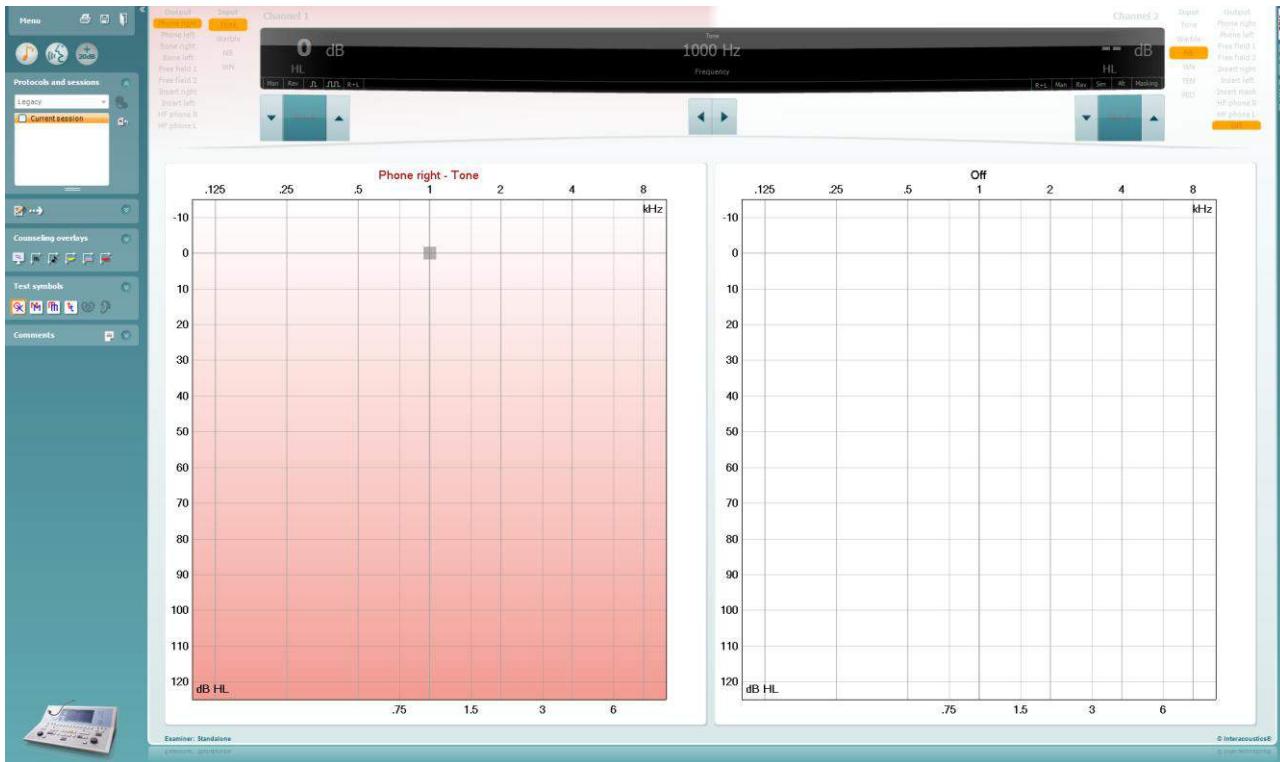
## 5.2

## 5.2 Transferring Measurement Data

When the instrument setup is done the Diagnostic Suite is ready for retrieving the audiogram data from the selected instrument. The Diagnostic Suite main display looks like this:

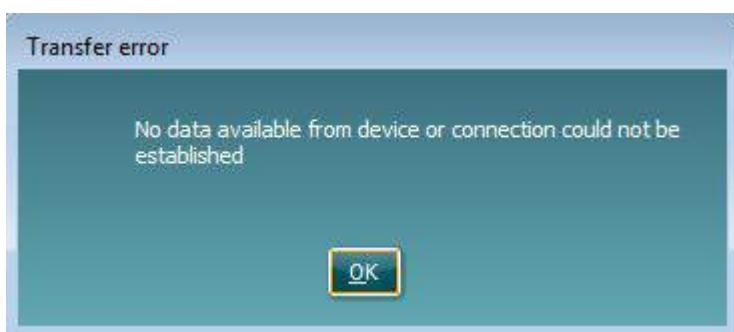
For the older IA instruments, the DS application version is for transferring data from the dedicated instruments only. The top display is disabled. For the newer hybrid audiometers (AD629/AC40/AD226), the instruments can be controlled by the suite. See next chapter.

When the measurement is completed on the selected instrument press the following icon to transfer the data:

**Important Note:** If the instrument is disconnected the following dialog is displayed when pressing the transfer button:

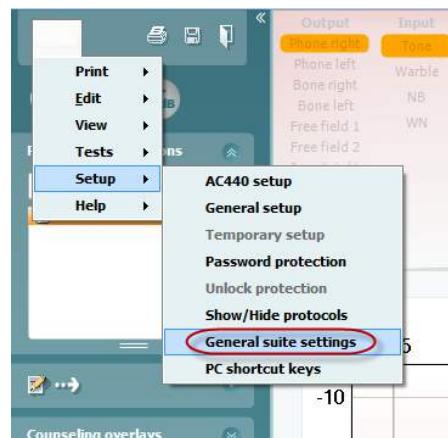
## 6 Diagnostic Suite with the New AD226/AD229/AD629/AC40/AT235(h)

This section describes the data-transfer and hybrid mode (On-Line / PC-operated modes) supported by the new AD226/AD229/AD629/AC40/AT235(h).

### 6.1 Instrument setup

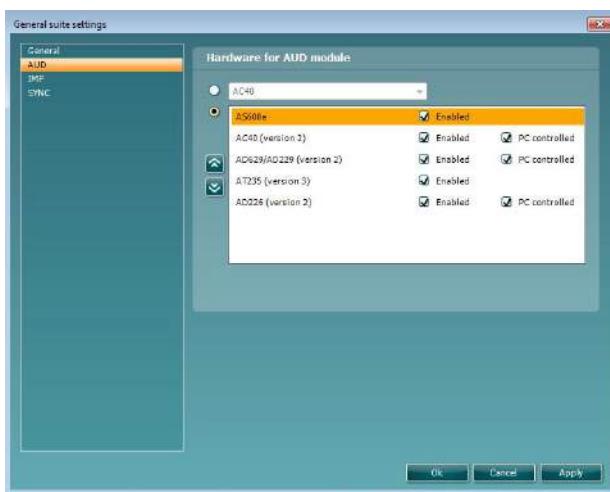
The instrument setup is done in the Diagnostic Suite under

**Menu | Setup | General suite settings** under the **AUD / IMP tab** section:

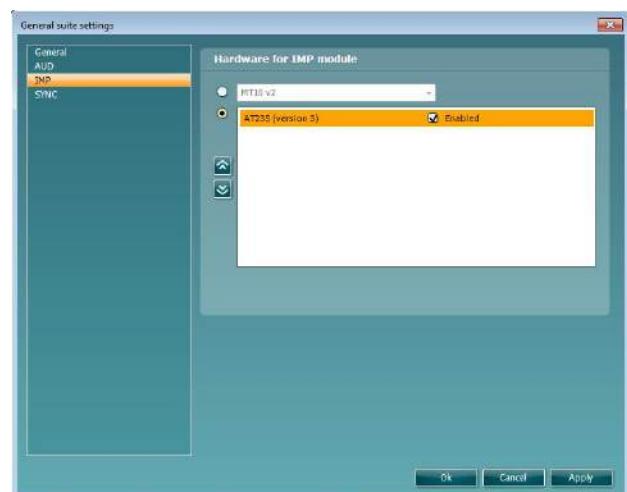
To setup AD226/AD229/AD629/AC40 select the AUD tab

To setup the New AT235 select IMP tab – please notices: the New AT235(h) can transfer tympanograms/audiograms and sync patient data. It's not possible use the AT235(h) in hybrid mode.

Hardware selection for AUD-module

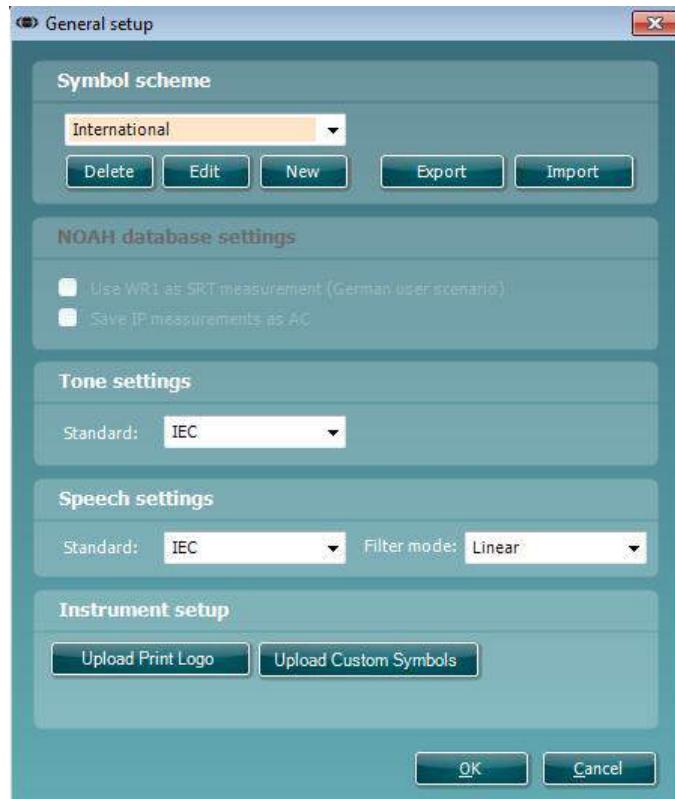
Hardware selection for IMP-module

**Important:** Please be sure to selected the “AD226 (version 2)”, “AD629 (version 2)”, “AC40 (version 2)” or “AT235 (version 2)” (and not the others, which refers to the old version).

**PC controlled instrument:** Unselect this if you want to run the AC40/AD629/AD226 as a standalone audiometer (i.e. not as a hybrid audiometer) but still being connected to the Diagnostic Suite. When pressing **Save Session** on the instrument, the session will automatically be transferred to the Diagnostic Suite. See below section “Sync Mode”.

## 6.2 General setup

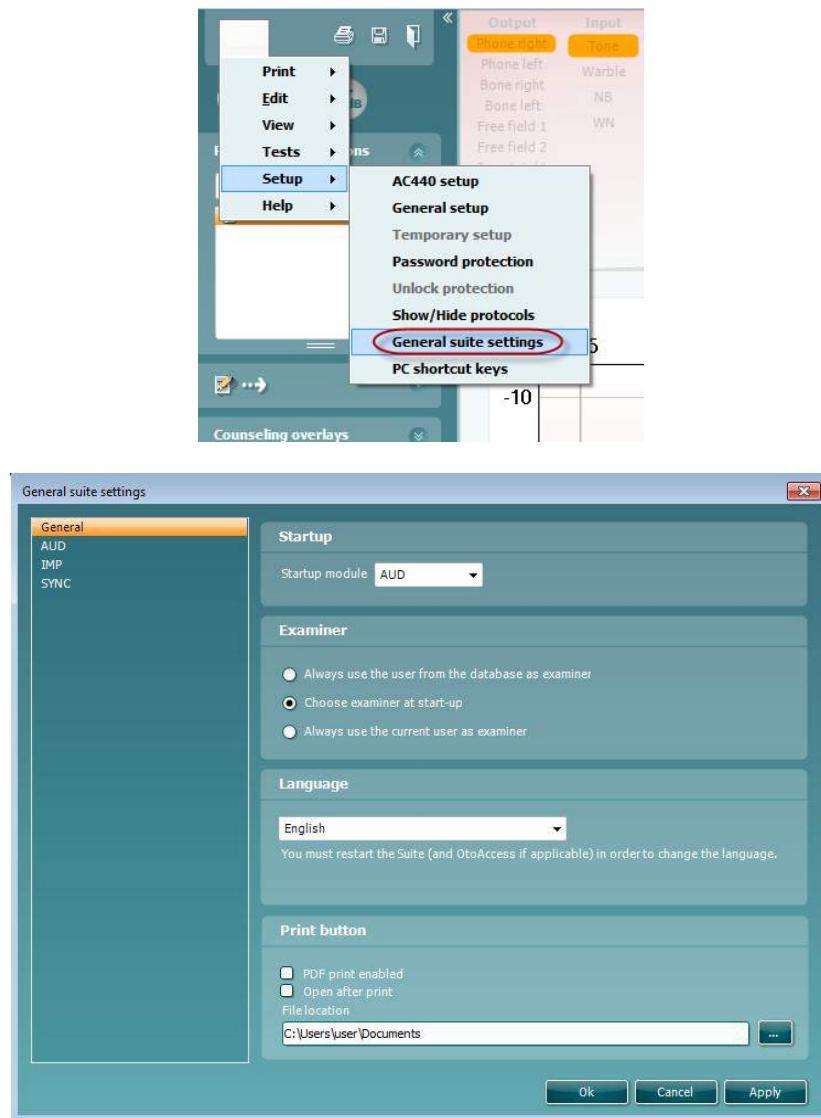
**Upload Print Logo and Audiogram Symbols:** A logo for direct print outs can be transferred to the New AC40/AD629/AT235(h) using the “Up Print Logo” button. The symbol scheme used in the Diagnostic Suite can be transferred to the AC40/AD629/AT235(h) (when viewing the build in audiogram) by using the “Upload Custom Symbols” button. Please refer to the AC40/AD629/AT235(h) operational manual for info on how to change the symbol scheme on the AC40/AD629/AT235(h).

**Changing standards:** after changing tone or speech standards the suite must be restarted.

## 6.3 Suite Settings

Menu | Setup | General suite settings under the AUD / IMP tab section:

**Startup:** you can select which module must be displayed when opening the suite

**Language:** select language. The suite must be restarted in order to change the language.

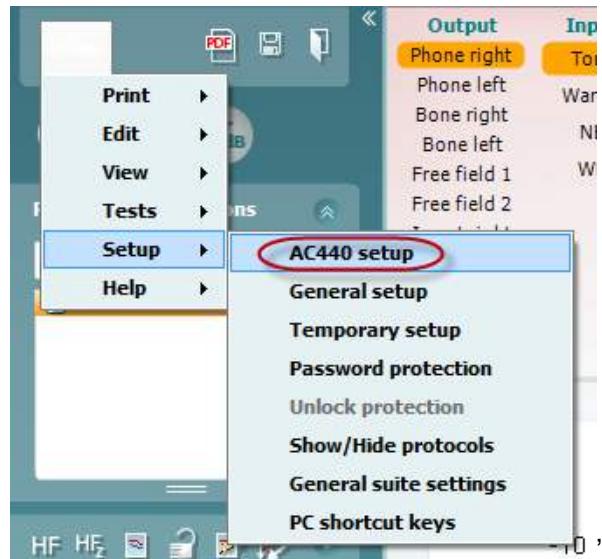
**Print button:** when “PDF printer enabled” is selected the print icon on the front screen will activate the PDF

printer.

## 6.4 Protocol setup

The protocol settings of the Diagnostic Suite AUD module can be modified under the AC440 setup:

Please refer to Additional Information document for more information about protocol configuration.

## 6.5 SYNC Mode

### 6.5.1 One click data transfer (Hybrid Mode disabled)

If the “PC controlled instrument” setting in the General Setup (see above) is deselected, the current audiogram will be transferred to the Diagnostic Suite as follows: When pressing *Save Session* on the instrument, the session will automatically be transferred to the Diagnostic Suite Start the suite with device connected

### 6.5.2 The Sync Tab

If there are several sessions stored on the instrument (under one or more patients) then the Sync tab must be used. The screen shot below shows the Diagnostic Suite with the SYNC tab open (underneath the AUD and IMP tabs in the upper right corner).

The SYNC tab provides the following possibilities:

**Client upload** is used for uploading clients from the database (Noah or OtoAccess™) to the audiometer.

**Session download** is used to download sessions (audiogram data) stored in the instrument memory into Noah, OtoAccess™ or XML (when running Diagnostic suite without a database).

### 6.5.3 Client Upload

The following screen shot shows the client upload screen:

- On the left side it's possible to search for the client in the database to transfer to the database using different search criteria's. Use the "Add" button to transfer (upload) the client from the database to the internal instrument memory.
- On the right side the clients currently stored on the internal instrument memory (hardware) is down. It's possible to remove all clients for individual clients using the "Remove all" or "Remove" buttons.

#### 6.5.4 Session download

The following screen shot shows the session download screen:

The screenshot shows the 'Session(s) on AD629 (Tone and Speech only)' screen. On the left, there's a sidebar with 'Client upload' and 'Session download' buttons. The main area displays a table with columns: Id, First name, Last name, Session(s), Status, and Action. There are five rows of data:

| Id  | First name | Last name | Session(s)   | Status          | Action        |
|-----|------------|-----------|--|-----------------|---------------|
|     | NoName     |           | 27. august 2012 14:53<br>27. august 2012 14:47<br>27. august 2012 14:45<br>27. august 2012 14:45<br>27. august 2012 14:44<br>27. august 2012 14:44<br>27. august 2012 14:43<br>27. august 2012 14:28 | No match (Skip) | <b>Change</b> |
| 7   | qq         | q         | 27. august 2012 14:47  | No match (Skip) | <b>Change</b> |
| 123 | Joan       | Jones     | 27. august 2012 14:46<br>2. august 2012 14:31  | No match (Skip) | <b>Change</b> |
| 777 |            |           | 22. august 2012 12:44<br>16. august 2012 13:51   | No match (Skip) | <b>Change</b> |

A large green 'Transfer to database' button is located at the top right of the table area. A small 'Info' icon is also present in the top right corner.

When pressing the icon the functionality of the "Session download" screen is described:

| Status                   | Meaning   |
|--------------------------|---|
| <b>Match (Transfer)</b>  | This client on AC40 (version 2) was found (matched) in the database and the measurement will be transferred (downloaded) into the database after pressing 'Transfer to database'.             |
| <b>No match (Skip)</b>   | This client on AC40 (version 2) was not found (not matched) in the database and the measurement will not be transferred (downloaded) into the database after pressing 'Transfer to database'. |
| <b>Download complete</b> | The client measurement data stored on AC40 (version 2) was successfully transferred (downloaded) to the selected client in the database.  |

A client on the AC40 (version 2) can be transferred (downloaded) into a different (existing or new) client in the database by selecting "Change" under the "Action" column. This will open a new dialog for changing the client selection.

## 6.6 Quick Start – Transferring and Saving Tympanometry Data

Diagnostic Suite allows you to transfer data, view data, add a report and save your data and report and/or print it with a customized print template.

To transfer tympanometry data you need to be in the **IMP** tab. At startup it will show empty graphs like below.

Press the arrow icon, in order to transfer data from your equipment to your pc.

On the New AT235(h) the transfer can also be done by pressing ‘Save session’ on the unit.

When the equipment is not yet connected or switched on or if the instrument setup is incorrect, a message will inform you to check the communication port.

Press the report icon, , in order to open the report editor and add notes to the transferred data.

Press the print icon, , in order to open the print wizard from which you can select according to which template to print.

Press the save icon, , to save or the save and exit icon, , to save and exit.

## 6.7 Hybrid (Online/PC-Controlled) Mode

### 6.7.1 Using the Tone Screen

The following section describes the elements of the tone screen.

**Menu** provides access to Print, Edit, View, Tests, Setup, and Help

**Print** allows for printing the session's acquired data.

**Save & New Session** saves the current session in Noah or OtoAccess™ and opens a new one.

**Save & Exit** saves the current session in Noah or OtoAccess™ and exits the Suite.

**Collapse** the left side panel.

**Go to Tone Audiometry** activates the tone screen when in another test.

**Go to Speech Audiometry** activates the speech screen when in another test.

**Extended Range +20 dB** extends the testing range and can be activated when the testing dial setting gets within 55 dB of the maximum level of the transducer.

Note that the extended range button will flash when it needs activation for reaching higher intensities.

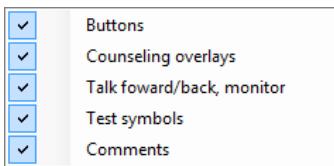
To switch on the extended range automatically, select the **Switch extended range on automatically** by going to the setup menu.

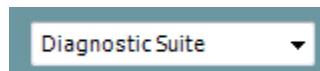
**Fold** an area so that it only shows the label or the buttons of that area.

**Unfold** an area so that all buttons and labels are visible

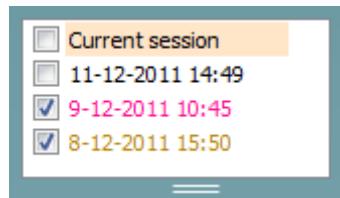
**Show/hide areas** can be found by right mouse clicking on one of the areas. The visibility of the different areas as well as the space that they take on the screen is locally saved to the examiner.

**List of Defined Protocols** allows for selecting a test protocol for the current test session. Right mouse click on a protocol allows the current examiner to set or deselect a default startup protocol. Please refer to “Additional Information” document for more information on protocols and protocol setup.

**Temporary Setup** allows for making temporary changes to the selected protocol. The changes will be valid for the current session only. After making the changes and returning to the main screen, the name of the protocol will be followed by an asterisk (\*).

**List of historical sessions** accesses historical sessions for comparison purposes. The audiogram of the selected session, indicated by the orange background, is shown in colours as defined by the used symbol set. All other audiograms that are selected by check marks show on screen in the colours as indicated by the text colour of the date and time stamp. Note that this listing can be resized by dragging the double lines up or down.

**Go to Current Session** brings you back to the current session.

**HF** High frequency

**High Frequency** shows frequencies on the audiogram (up to 20 kHz for the AC40/AD629). However, you will only be able to test in the frequency range the selected headset is calibrated for.

**HF<sub>z</sub>** High frequency zoom

**High Frequency Zoom**<sup>1</sup> activates high frequency testing and zooms in on the high frequency range.

**Single audiogram**

**Single audiogram** toggles between viewing the information of both ears in a single graph and two separate graphs.

**MF** Multi frequencies

**Multi frequencies**<sup>2</sup> activate testing with frequencies in between the standard audiogram points. The frequency resolution can be adjusted in the AC440 setup.

**Synchronize channels**

**Synchronize channels** locks the two channels together. This function may be used to perform synchronous masking.

**Edit mode**

**Edit Mode** button activates the editing function. Left clicking on the graph will add/move a point to the position of the cursor. If right clicking on a specific stored point a context menu appears offering you the following options:

- Add unmasked threshold
- Add no response
- Add masked threshold
- Add masked-no-response threshold
- Copy bone thresholds to other ear
- Delete threshold
- Delete curve
- Hide unmasked thresholds where masked exist
- Cancel

Intensity: 15dB Frequency: 250Hz

**Mouse controlled audiology**

**Mouse controlled audiology** enables you to do the audiology using the mouse only. Left click on the mouse to present the stimulus. Right click on the mouse to store the result.

**dB step size**

The **dB step size** button indicates to which dB step size the system is currently set. It rotates between 1 dB, 2 dB and 5 dB step sizes.

**Hide unmasked thresholds**

The **hide unmasked threshold** will hide those unmasked thresholds where masked thresholds exist.

**Toggle masking help**

**Toggle Masking Help** will activate or deactivate the Masking Help feature.

For more information on Masking Help, please refer to the “Additional Information” or the “Masking Help Quick Guide” documents.

<sup>1</sup> HF requires an additional license for the AC440. If not purchased, the button is grayed out.

<sup>2</sup> MF requires additional license for the AC440. If not purchased, the button is grayed out.

**Toggle Automasking** will activate or deactivate the Automasking feature. For more information on Automasking, please refer to the “Additional Information” or the “Masking Help Quick Guide” documents

**Talk Forward** activates the Talk Forward microphone. The arrow keys can be used to set the talk forward level through the currently selected transducers. The level will be accurate when VU meter indicates to be at zero dB.

Selecting the **Monitor Ch1** and/or **Ch2** check boxes allow you to monitor one or both channels through an external loudspeaker/headset connected to the monitor input. The monitor intensity is adjusted by the arrow keys.

The **Talk back** check box enables you to listen to the patient. Note that you need to be equipped with a microphone connected to the talk back input and an external loudspeaker/headset connected to the monitor input.

The **Patient monitor** opens an always-on-top window with the tone audiograms and all its counseling overlays shown. The size and position of the patient monitor gets saved for each examiner individually.

The **Phonemes** counseling overlay shows phonemes as it is set up in the protocol that is currently in use.

The **Sound examples** counseling overlay shows pictures (png-files) as they are set up in the protocol that is currently in use.

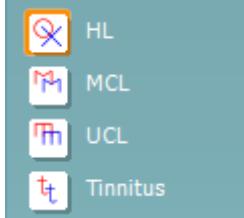
The **Speech banana** counseling overlay shows the speech area as it is set up in the protocol that is currently in use.

The **Severity** counseling overlay shows the degrees of hearing loss as it is set up in the protocol that is currently in use.

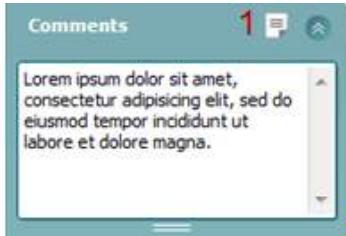
The **Max. testable values** shows the area beyond the maximum intensity the system allows. This is a reflection of the transducer calibration and depends on the extended range being activated.

Selecting **HL**, **MCL**, **UCL** or **Tinnitus** sets the symbol types that are currently in use by the audiogram. HL stands for hearing level, MCL stands for most comfortable level and UCL stands for uncomfortable level. Note that these buttons show the unmasked right and left symbols of the currently selected symbol set.

Each type of measurement is saved as a separate curve.

In the **Comments** section you can type comments related to any audiometric test. The used space by the comments area can be set by dragging the double line with your mouse. Pressing the **Report editor** (1) button opens a separate window for adding notes to the current session. The report editor and comment box contain the same text. In case the formatting of the text is important, this can only be set within the report editor.

After saving the session, changes can only be made within the same day until the date changes (at midnight). **Note:** these timeframes are limited by HIMA and the Noah software, and not by Interacoustics.

| Output       | Input  |
|--------------|--------|
| Phone right  | Tone   |
| Phone left   | Warble |
| Bone right   | NB     |
| Bone left    | WN     |
| Free field 1 |        |
| Free field 2 |        |
| Insert right |        |
| Insert left  |        |

The **Output** list for channel 1 provides the option to test through head phones, bone conductor, free field speakers or insert phones. Note that the system only shows the calibrated transducers.

The **Input** list for channel 1 provides the option to select pure tone, warble tone, narrow band noise (NB) and white noise (WN).

Note that the background shading is according to the side that is selected, red for right and blue for left.

| Input  | Output       |
|--------|--------------|
| Tone   | Phone right  |
| Warble | Phone left   |
| NB     | Free field 1 |
| WN     | Free field 2 |
| TEN    | Insert right |
|        | Insert left  |
|        | Insert mask  |
|        | Off          |

The **Output** list for channel 2 provides the option to test through head phones, free field speakers, insert phones or insert masking phone. Note that the system only shows the calibrated transducers.

The **Input** list for channel 2 provides the option to select pure tone, warble tone, narrow band noise (NB), white noise (WN) and TEN noise<sup>3</sup>.

Note that the background shading is according to the side that is selected, red for right, blue for left, and white when off.

**Pulsation** allows for single and continuous pulsating presentation. The duration of the stimulus can be adjusted in the AC440 setup.

**Sim/Alt** allows toggling between **Simultaneous** and **Alternate** presentation. Ch1 and Ch2 will present the stimulus simultaneously when Sim is selected. When Alt is selected, the stimulus will alternate between Ch1 and Ch2.

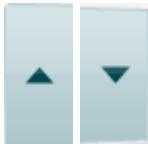
**Masking** indicates if channel 2 is currently in use as a masking channel and in that way makes sure masking symbols are used in the audiogram. For example in paediatric testing through free field speakers, channel 2 can be set as a second testing channel. Note that a separate store function for channel 2 is available when channel 2 is not used for masking.

The **Right + Left** allows for presenting tones in both ears in channel 1 and noise in both ears in channel 2.

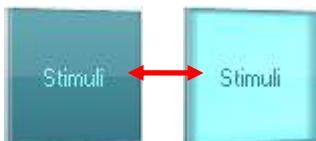
<sup>3</sup> TENS test requires an additional license for the AC440. If not purchased, the stimulus is grayed out.

**dB HL Increase** and **Decrease** buttons allows for increasing and decreasing the intensities of channel 1 and 2.

The arrow keys on the PC keyboard can be used for in-/decreasing channel 1 intensities.

PgUp and PgDn on the PC keyboard can be used for in-/decreasing channel 2 intensities.

**Stimuli** or **attenuator** buttons will light up when the mouse goes over and indicates the active presentation of a stimulus.

A right mouse click in the Stimuli area will store a no response threshold. A left mouse click in the Stimuli area will store the threshold at the current position.

Channel 1 stimulation can also be obtained by pressing the space bar or left Ctrl key on the PC keyboard.

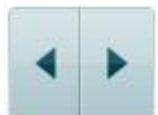
Channel 2 stimulation can also be obtained by pressing the right Ctrl key on the PC keyboard.

Mouse movements in the Stimuli area for both channel 1 and channel 2 can ignored depending on the setup.

**Frequency and Intensity display** area shows what is currently presented. To the left the dB HL value for channel 1 is shown and to the right for channel 2 In the centre the frequency is displayed.

Notice that the dB dial setting will flash when trying to go louder than the maximum available intensity.

**Frequency increase/decrease** increases and decreases the frequency respectively. This can also be obtained using the left and right arrow keys on the PC keyboard.

**Storing** thresholds for channel 1 is done by pressing **S** or by a left mouse click in the Stimuli button of channel 1. Storing a no response threshold can be done by pressing **N** or by a right mouse click on the Stimuli button of channel 1.

**Storing** thresholds for channel 2 is available when channel 2 is not the masking channel. It is done by pressing **<Shift> S** or by a left mouse click on the Stimuli button of channel 2. Storing a no response threshold can be done by pressing **<Shift> N** or by a right mouse click in the attenuator of channel 2.

**The hardware indication picture** indicates whether the hardware is connected. **Simulation mode** is indicated when operating the software without hardware.

When opening the Suite the system will automatically search for the hardware.

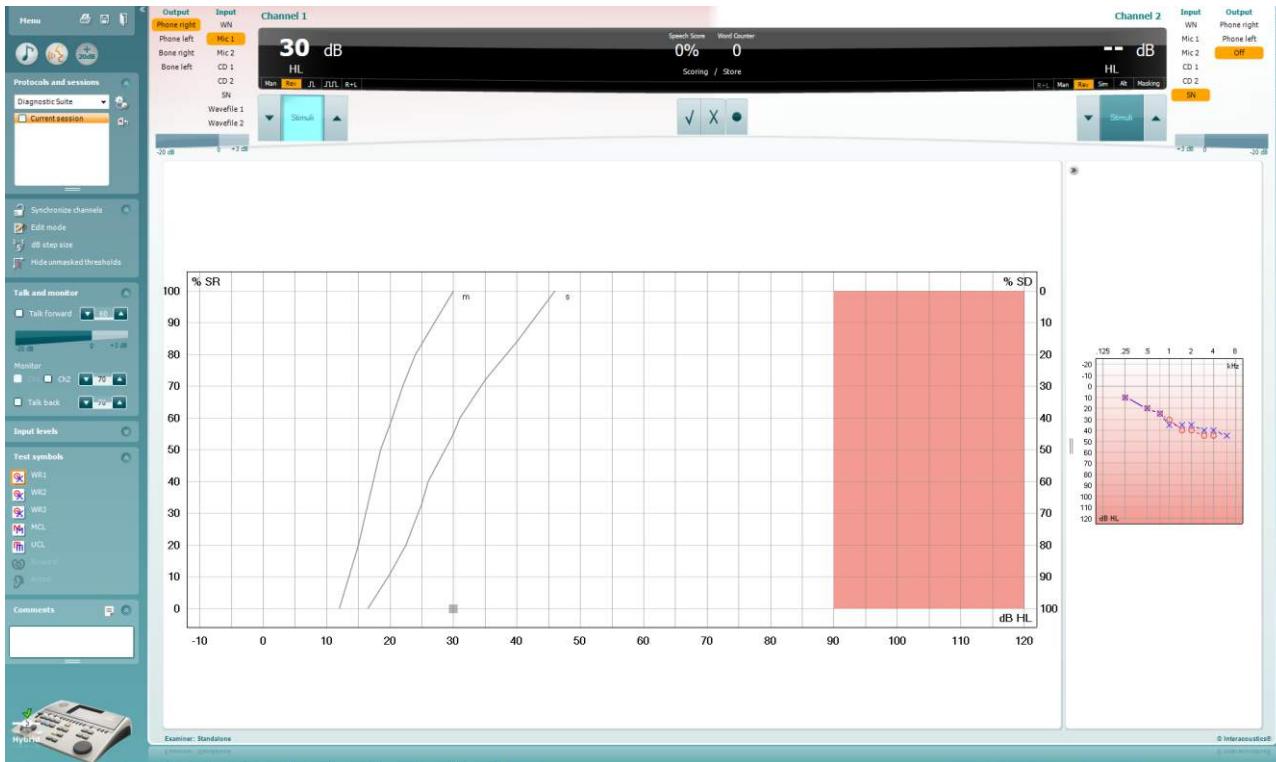
The **Examiner** indicates the current clinician who is testing the patient. The examiner is saved with a session and can be printed with the results.

For each examiner is logged how the suite is set up with regards to the use of space in the screen. The examiner will find that the suite starts up looking the same as the last time they used the software. An examiner can also select which protocol must be selected at start up (by right mouse click on the protocol selection list).

### 6.7.2 Using the Speech Screen

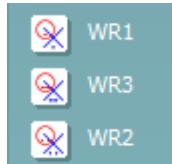
The following section describes the elements of the speech screen in addition to the tone screen:

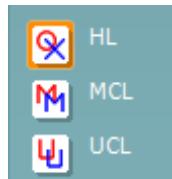
**Input levels** sliders allow for adjusting the input level to 0 VU for the selected input. This ensures that correct calibration is obtained for Mic1, Mic2, CD1, and CD2.

|      | ▲  |
|------|----|
| Mic1 | 28 |
| Mic2 | 27 |
| CD1  | 26 |
| CD2  | 26 |

**WR1, WR2 and WR3 (Word Recognition)** allows selecting different speech list setups as defined by the selected protocol. The labels of these lists which go along with these buttons can also be customized in the protocol setup.

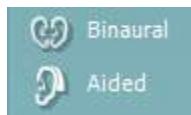
Selecting **HL, MCL and UCL** sets the symbol types that are currently in use by the audiogram. HL stands for hearing level, MCL stands for most comfortable level and UCL stands for uncomfortable level.

Each type of measurement is saved as a separate curve.

**Binaural** and **Aided** function allows for indicating if the test is performed binaurally or while the patient is wearing hearing aids. This feature is only active in the Speech Audiometry screen.

| Output       | Input      |
|--------------|------------|
| Phone right  | WN         |
| Phone left   | Mic 1      |
| Bone right   | Mic 2      |
| Bone left    | Cd 1       |
| Free field 1 | Cd 2       |
| Free field 2 | SN         |
| Insert right | Wavefile 1 |
| Insert left  | Wavefile 2 |

The **Output** list for channel 1 provides the option to test through the desired transducers. Note that the system only shows the calibrated transducers.

The **Input** list for channel 1 provides the option to select white noise (WN), speech noise (SN), microphone 1 or 2 (Mic1 and Mic2), CD1, CD2 and wavefile.

Note that the background shading is according to the side that is selected, red for right and blue for left.

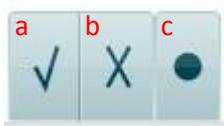
The **Output** list for channel 1 provides the option to test through the desired transducers. Note that the system only shows the calibrated transducers.

| Input      | Output       |
|------------|--------------|
| WN         | Phone right  |
| Mic 1      | Phone left   |
| Mic 2      | Free field 1 |
| Cd 1       | Free field 2 |
| Cd 2       | Insert right |
| SN         | Insert left  |
| Wavefile 1 | Insert mask  |
| Wavefile 2 | Off          |

The **Input** list for channel 2 provides the option to select white noise (WN), speech noise (SN), microphone (Mic1 and Mic2), CD1, CD2 and wavefile.

Note that the background shading is according to the side that is selected, red for right, blue for left, and white when off.

### Speech Scoring:

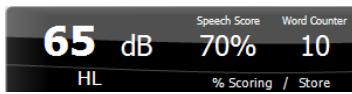
- a) **Correct:** A mouse click on this button will store the word as correctly repeated. You can also click on the **Left** arrow key to store as correct.
- b) **Incorrect:** A mouse click on this button will store the word as incorrectly repeated. You can also click on the **Right** arrow key to store as incorrect
- c) **Store:** A mouse click on this button will **store** the speech threshold in the speech graph. A point can also be stored by pressing **S**.

### Phoneme scoring:

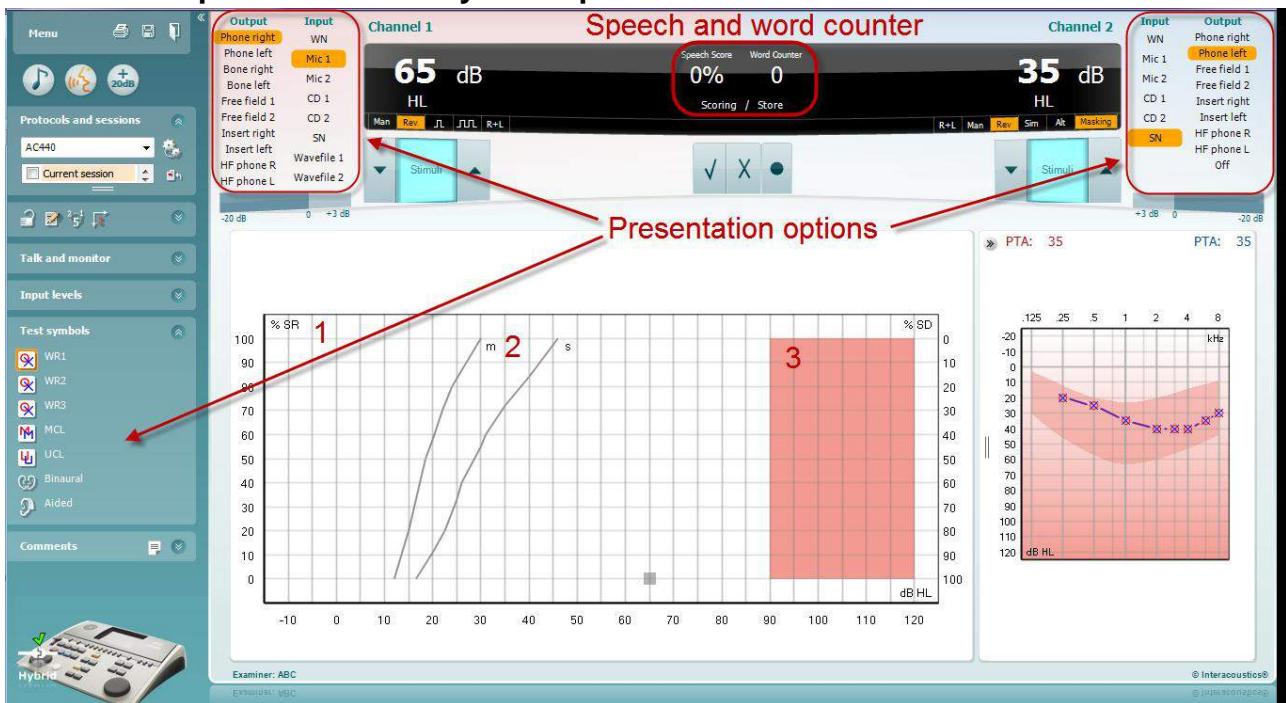
- a) **Phoneme scoring:** If phoneme scoring is selected in the AC440 setup, mouse click on the corresponding number to indicate phoneme score. You can also click on the **Up** key to store as correct and **Down** key to store as incorrect.
- b) **Store:** A mouse click on this button will store the speech threshold in the speech graph. A point can also be stored by pressing **S**.

**Frequency and Speech score display** shows what is currently presented. On the left the dB value for channel 1 is shown and on the right side for channel 2.

In the centre of the current *Speech Score* in % and the *Word Counter* monitors the number of words presented during the test.

### 6.7.3 Speech Audiometry in Graph Mode

Graph mode presentation settings under “Test Symbols” and in the presentation options (Ch1 and Ch2) in the upper part of the screen shows where you can adjust the test parameters during the test.

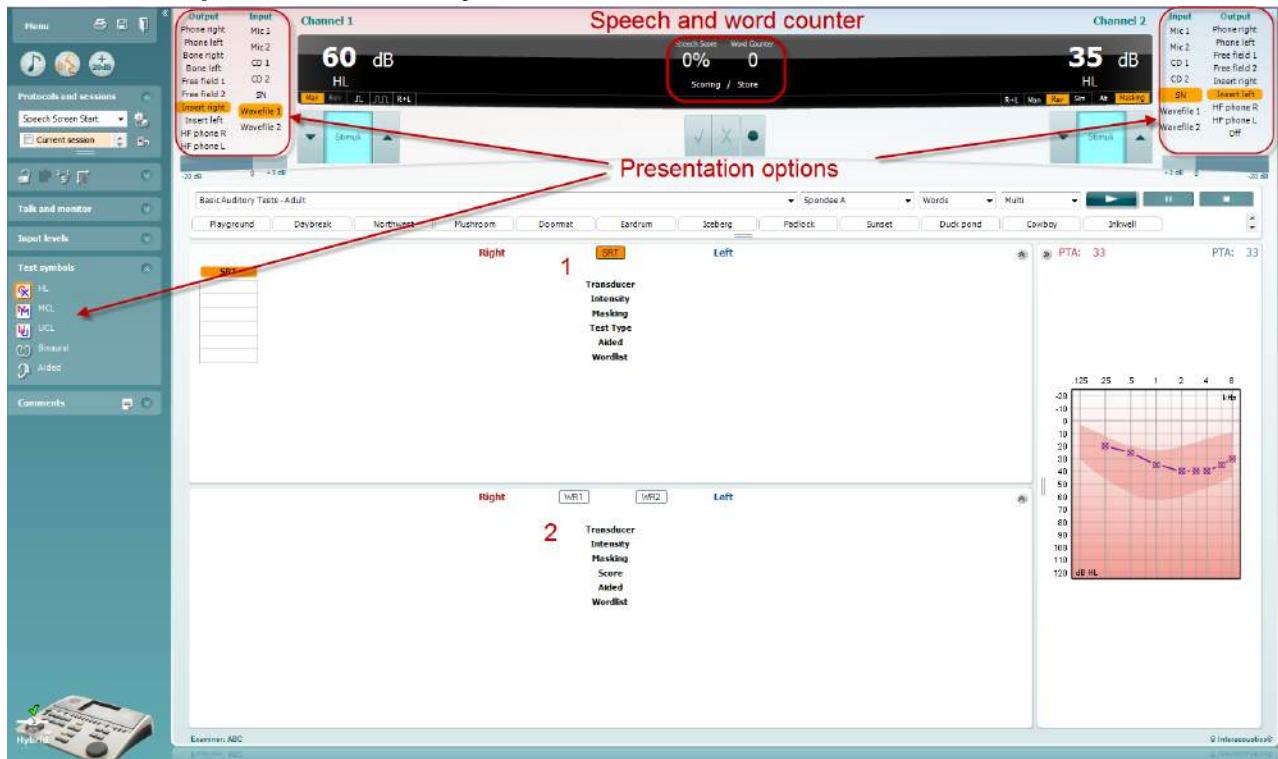
- 1) **The graph:** The curves of the recorded speech graph will be displayed on your screen.

The x-axis shows the intensity of the speech signal and the y-axis shows the score in percent.

The score is also displayed in the black display in the upper part of the screen, along with a word counter.

- 2) **The norm curves** illustrate norm values for **S** (Single syllabic) and **M** (Multi syllabic) speech material respectively. The curves can be edited according to individual preferences in the AC440 setup .
- 3) **The shaded area** illustrates the maximum intensity the system will allow. The *Extended Range +20 dB* button can be pressed to go higher. The maximal loudness is determined by the transducer calibration.

### 6.7.4 Speech Audiometry in Table Mode

The AC440 Table Mode consists of two tables:

- 1) The **SRT** (Speech Reception Threshold) table. When the SRT test is active, it is indicated in orange **SRT**
- 2) The **WR** (Word Recognition) table. When WR1, WR2, or WR3 is active the corresponding label will be orange **WR1**

### The SRT table

The SRT table (Speech Reception Threshold table) allows for measuring multiple SRTs using different test parameters, e.g. *Transducer*, *Test Type*, *Intensity*, *Masking*, and *Aided*.

Upon changing *Transducer*, *Masking*, and/or *Aided* and re-testing, an additional SRT entry will appear in the SRT table. This allows for multiple SRT measurements to be shown in the SRT table.

Please refer to [Additional Information](#) document for more information about SRT testing.

| Right     |           | SRT        | Left      |           |
|-----------|-----------|------------|-----------|-----------|
| SRT       | SRT       | Transducer | SRT       | SRT       |
| Phone     | Phone     | Intensity  | Phone     | Phone     |
| 30        | 10        | Masking    | 10        | 30        |
| 15        | 15        | Test Type  | 15        | 15        |
| HL        | HL        | Aided      | HL        | HL        |
| x         | x         | Wordlist   | x         |           |
| Spondee A | Spondee B |            | Spondee A | Spondee B |

### The WR Table

The word recognition (WR) table allows for measuring multiple WR scores using different parameters (e.g. *Transducer*, *Test Type*, *Intensity*, *Masking*, and *Aided*).

Upon changing Transducer, Masking, and/or Aided re-testing an additional WR entry will appear in the WR table. This allows for multiple WR measurements to be shown in the WR table.

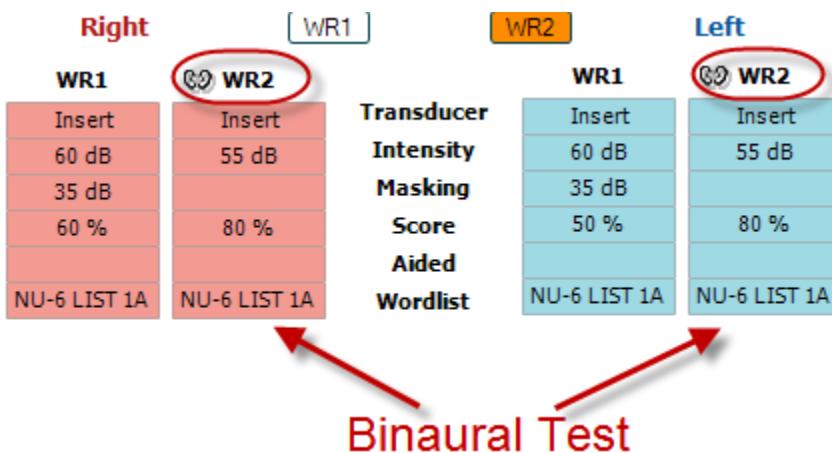
Please refer to the Callisto Additional Information document for more information about Word Recognition testing.

| Right        |              | WR1        | WR2          | WR3       | Left |
|--------------|--------------|------------|--------------|-----------|------|
| WR1          | WR1          | Transducer | WR1          | WR2       |      |
| Phone        | FF1          | Intensity  | Phone        | FF2       |      |
| 55           | 55           | Masking    | 55           | 30        |      |
| 85           | 95           | Score      | 90           |           |      |
| x            | x            | Aided      | 100          |           |      |
| NU-6 LIST 1A | NU-6 LIST 3A | Wordlist   | NU-6 LIST 1A | Spontee A |      |

### Binaural and Aided options

To perform binaural speech tests:

1. Click on either SRT or WR, to choose the test to be conducted binaurally
2. Ensure that the transducers are set up for binaural testing. For example, insert Right in channel 1 and insert Left in channel 2
3. Click on  Binaural
4. Proceed with the test; when stored, results will be stored as binaural results

| Right        | WR1          | WR2          | Left         |
|--------------|--------------|--------------|--------------|
| WR1          | WR2          | WR2          | WR2          |
| Insert       | Insert       | Insert       | Insert       |
| 60 dB        | 55 dB        | 60 dB        | 55 dB        |
| 35 dB        |              | 35 dB        |              |
| 60 %         | 80 %         | 50 %         | 80 %         |
| NU-6 LIST 1A | NU-6 LIST 1A | NU-6 LIST 1A | NU-6 LIST 1A |

**Binaural Test**

To perform an aided test:

1. Select the desired transducer. Typically, aided testing is done in the Free Field. However, in certain conditions, it could be possible to test deeply inserted CIC hearing instruments under headphones, which would show ear-specific results
2. Click on the Aided button
3. Click on the Binaural button if the test is done in the Free Field so that the results are stored for both ears at the same time
4. Proceed with test; results will then be stored as aided by showing an Aided icon

| WR2   |
|---|
| FF1   |
| 15 dB   |
| 80 %  |
|  |
| NU-6 LIST 3A  |

### 6.7.5 PC Keyboard Shortcuts Manager

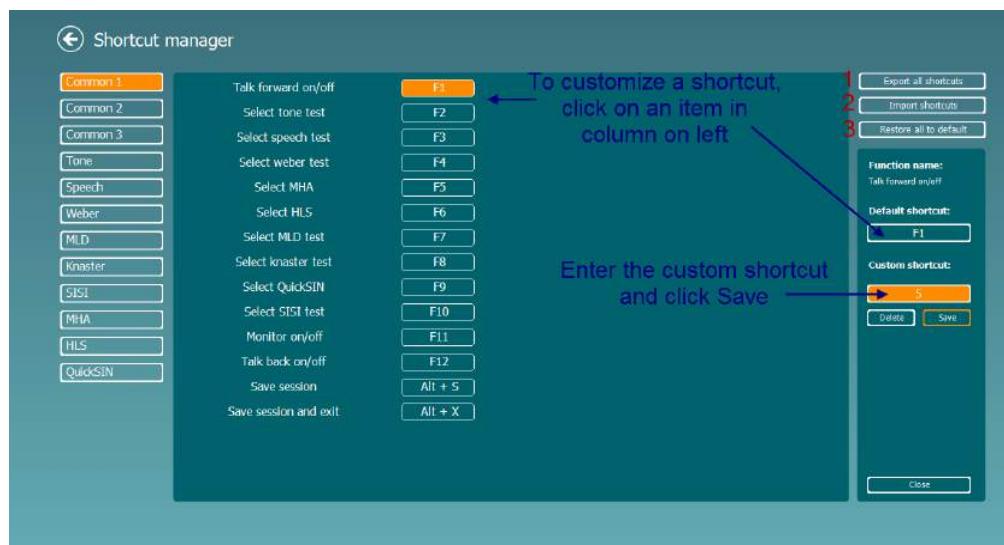
The PC Shortcut Manager allows the user to personalize PC shortcuts in the AUD Module. To access the PC Shortcut Manager:

**Go to AUD module | Menu | Setup | PC Shortcut Keys**

To view the default shortcuts, click on the items in the left hand column (Common 1, Common 2, Common 3, etc.)

To personalize a shortcut, click on the column in the middle and add the custom shortcut in the field on the right hand of the screen

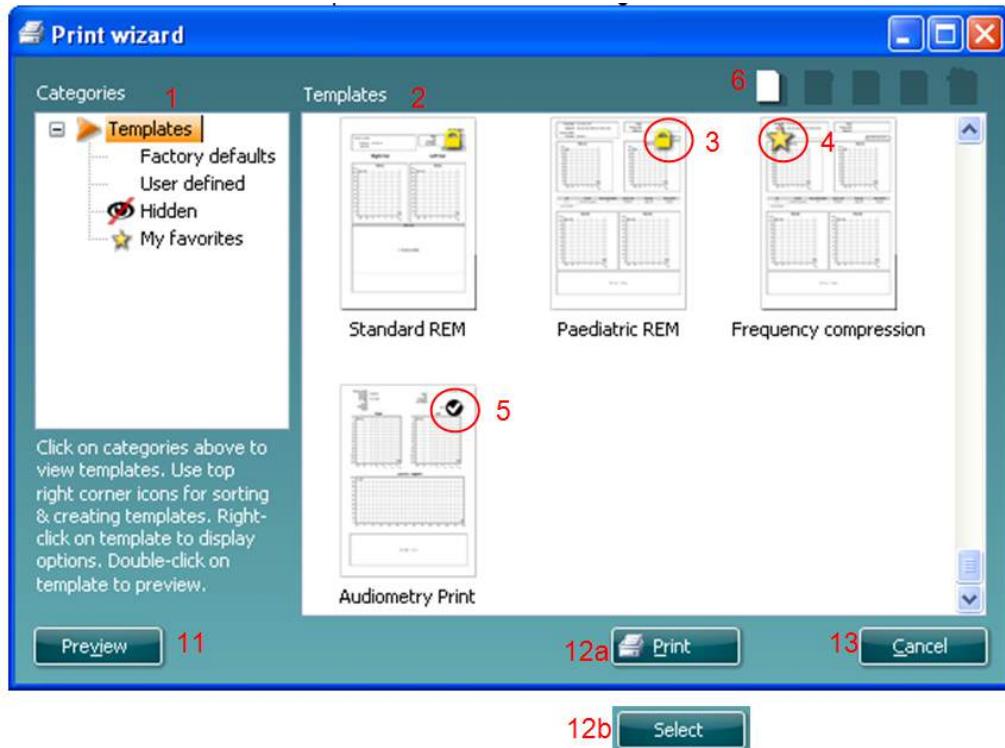
1. **Export all shortcuts:** Use this function to save custom shortcuts and transfer them to another computer
2. **Import shortcuts:** Use this function to import shortcuts that have already been exported from another computer
3. **Restore all defaults:** Use this function to restore the PC shortcuts to Factory Settings default

## 7 Using the Print Wizard

In the Print Wizard you have the option to create customized print templates which can be linked to individual protocols for quick printing. The Print Wizard can be reached in two ways.

- If you want to use a template for general use, or select an existing one for printing: Go to **Menu/ File/Print Layout...** in both AUD and IMP tabs.
- If you want to create a template or select an existing one to link to a specific AUD protocol: Select a specific protocol and select **Menu | Setup | AC440 setup**. Select the specific protocol from the drop down menu and select **Print Setup** at the bottom of the window.

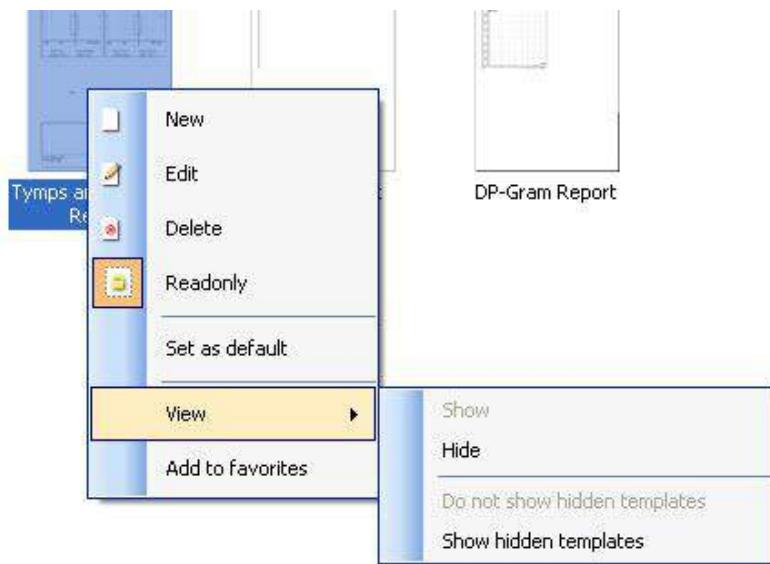
Now the **Print Wizard** window opens and shows the following information and functionalities:

- Underneath **Categories** you can select
  - Templates** to show all available templates
  - Factory defaults** to show only standard templates
  - User defined** to show only custom templates
  - Hidden** to show hidden templates
  - My favorites** to show only templates marked as a favorite
- Available templates from the selected category are shown in the **Templates** viewing area.
- Factory default templates are recognized by the lock icon. They ensure that you always have a standard template and do not need to create a customized one. However, in order to edit these default templates, they need to be saved under a new name. **User defined**/created templates can be set to **Read-only** (showing the lock icon), by right clicking on the template and selecting **Read-only** from the drop down list. **Read-only** status can also be removed from **User defined** templates by following the same steps.
- Templates added to **My favorites** are marked with a star. Adding templates to **My favorites** allows quick viewing of your most commonly used templates.
- The template that is attached to the selected protocol when entering the print wizard is recognized by a checkmark.
- Press the **New Template** button to open a new empty template.

7. Select one of the existing templates and press the **Edit Template** button to modify the selected layout.
8. Select one of the existing templates and press the **Delete Template** button to delete the selected template. You will be prompted to confirm that you want to delete the template.
9. Select one of the existing templates and press the **Hide Template** button to hide the selected template. The template will now be visible only when **Hidden** is selected under **Categories**. To unhide the template, select **Hidden** under **Categories**, right click on the desired template and select **View>Show**.
10. Select one of the existing templates and press the **My Favorites** button to mark the template as a favorite. The template can now be quickly found when **My Favorites** is selected under **Categories**. To remove a template marked with a star from My Favorites, select the template and press the **My Favorites** button.
11. Select one of the templates and press the **Preview** button to print preview the template on screen.
12. Depending how you reached the Print Wizard, you will have the option to press
  - a. **Print** for using the selected template for printing or press
  - b. **Select** for dedicating the selected template to the protocol from which you got into the Print Wizard.
13. To leave the Print Wizard without selecting or changing a template press **Cancel**.

Right clicking on a specific template provides a drop down menu offering an alternative method for performing the options as described above:

For more information related to the Print reports and Print Wizard, please refer to the Additional Information document or the Print Report Quick Guide on [www.interacoustics.com](http://www.interacoustics.com)

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| Status  | Meaning   |
|---|---|
|  <b>Match (Transfer)</b> | This client on AC40 (version 2) was found (matched) in the database and the measurement will be transferred (downloaded) into the database after pressing 'Transfer to database'.             |
| <b>No match (Skip)</b>  | This client on AC40 (version 2) was not found (not matched) in the database and the measurement will not be transferred (downloaded) into the database after pressing 'Transfer to database'. |
| <b>Download complete</b>  | The client measurement data stored on AC40 (version 2) was successfully transferred (downloaded) to the selected client in the database.  |

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The screenshot displays three panels. The left panel is a sidebar with "Menu", "Client upload", and "Session download". The middle panel is titled "Client Transfer to AD629" and contains a "Client Search and Select" table with columns: Last name, First name, Birthdate, Id, and Address. It shows one row: "Standalone" under "Last name", "01-01-2008" under "Birthdate", and "01-01-2008" under "Address". The right panel is titled "Clients on hardware" and shows a table with columns: Last name, First name, and Id. It lists three clients: "Jones" (Id 123, First name Joan), "qq" (Id 777, First name NoName), and "q" (Id 7, First name NoName). A green arrow points from the "Client upload" button in the sidebar to the "Add" button in the search table.

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The screenshot shows a sidebar with "Menu", "Client upload", and "Session download". The main panel is titled "Session(s) on AD629 (Tone and Speech only)" and contains a table with columns: Id, First name, Last name, Session(s), Status, and Action. It lists four sessions:
 

| Id  | First name | Last name | Session(s)   | Status          | Action                  |
|-----|------------|-----------|--|-----------------|-------------------------|
|     |            |           |  |                 | Change                  |
|     | NoName     |           | 27. august 2012 14:53<br>27. august 2012 14:47<br>27. august 2012 14:45<br>27. august 2012 14:45<br>27. august 2012 14:44<br>27. august 2012 14:44<br>27. august 2012 14:43<br>27. august 2012 14:28 | No match (Skip) | <button>Change</button> |
| 7   | qq         | q         | 27. august 2012 14:47  | No match (Skip) | <button>Change</button> |
| 123 | Joan       | Jones     | 27. august 2012 14:46<br>2. august 2012 14:31  | No match (Skip) | <button>Change</button> |
| 777 |            |           | 22. august 2012 12:44<br>16. august 2012 13:51   | No match (Skip) | <button>Change</button> |

 A green arrow points from the "Session download" button in the sidebar to the "Transfer to database" button in the header of the main panel.

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| Status            | Meaning   |
|-------------------|---|
| Match (Transfer)  | This client on AC40 (version 2) was found (matched) in the database and the measurement will be transferred (downloaded) into the database after pressing 'Transfer to database'.             |
| No match (Skip)   | This client on AC40 (version 2) was not found (not matched) in the database and the measurement will not be transferred (downloaded) into the database after pressing 'Transfer to database'. |
| Download complete | The client measurement data stored on AC40 (version 2) was successfully transferred (downloaded) to the selected client in the database.  |

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The first screenshot shows the main menu with options for Client upload and Session download. A green arrow points to the Client upload icon. The second screenshot shows a search interface titled "Client Transfer to AD629" with a table containing one row: "Standalone" with ID "01-01-2008". The third screenshot shows a list titled "Clients on hardware" with three entries: Jones (NoName, Joan), qq, and 777.

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The screenshot shows the main menu with Client upload and Session download options. The Session download icon has a green arrow pointing to it. The main window displays a table titled "Session(s) on AD629 (Tone and Speech only)" with four rows of data. The columns are Id, First name, Last name, Session(s), Status, and Action. The data is as follows:

| Id  | First name | Last name | Session(s)   | Status          | Action |
|-----|------------|-----------|--|-----------------|--------|
|     | NoName     |           | 27. august 2012 14:53<br>27. august 2012 14:47<br>27. august 2012 14:45<br>27. august 2012 14:45<br>27. august 2012 14:44<br>27. august 2012 14:44<br>27. august 2012 14:43<br>27. august 2012 14:28 | No match (Skip) | Change |
| 7   | qq         | q         | 27. august 2012 14:47  | No match (Skip) | Change |
| 123 | Joan       | Jones     | 27. august 2012 14:46<br>2. august 2012 14:31  | No match (Skip) | Change |
| 777 |            |           | 22. august 2012 12:44<br>16. august 2012 13:51   | No match (Skip) | Change |

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The screenshot shows the 'Client Transfer to AD629' interface. On the left, there's a sidebar with icons for 'Client upload' and 'Session download'. The main area has two sections: 'Client Search and Select' containing a table with columns Last name, First name, Birthdate, Id, and Address, and 'Clients on hardware' containing a table with columns Last name, First name, and Id.

| Last name    | First name | Id  |
|--------------|------------|-----|
| <b>Jones</b> | NoName     | 123 |
| qq           | Joan       | 777 |

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The screenshot shows the 'Session(s) on AD629 (Tone and Speech only)' interface. It features a table with columns Id, First name, Last name, Session(s), Status, and Action. Below the table, there's a detailed view for a specific session.

| Id     | First name | Last name | Session(s)   | Status          | Action |
|--------|------------|-----------|--|-----------------|--------|
| NoName |            |           | 27. august 2012 14:53<br>27. august 2012 14:47<br>27. august 2012 14:45<br>27. august 2012 14:45<br>27. august 2012 14:44<br>27. august 2012 14:44<br>27. august 2012 14:43<br>27. august 2012 14:28 | No match (Skip) | Change |
| 7      | qq         | q         | 27. august 2012 14:47  | No match (Skip) | Change |
| 123    | Joan       | Jones     | 27. august 2012 14:46<br>2. august 2012 14:31  | No match (Skip) | Change |
| 777    |            |           | 22. august 2012 12:44<br>16. august 2012 13:51   | No match (Skip) | Change |

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| Status                   | Meaning   |
|--------------------------|---|
| Match (Transfer)         | This client on AC40 (version 2) was found (matched) in the database and the measurement will be transferred (downloaded) into the database after pressing 'Transfer to database'.             |
| <b>No match (Skip)</b>   | This client on AC40 (version 2) was not found (not matched) in the database and the measurement will not be transferred (downloaded) into the database after pressing 'Transfer to database'. |
| <b>Download complete</b> | The client measurement data stored on AC40 (version 2) was successfully transferred (downloaded) to the selected client in the database.  |

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